

Safety Data Sheet

Benzene

Division of Safety
National Institutes
of Health



WARNING!

THIS COMPOUND IS ACUTELY TOXIC, CARCINOGENIC, EMBRYOTOXIC, AND MUTAGENIC. IT IS READILY ABSORBED BY VARIOUS BODY TISSUES THROUGH THE SKIN AND RESPIRATORY AND INTESTINAL TRACTS AND TRANSPLACENTALLY. IT MAY IRRITATE THE EYES, MUCOUS MEMBRANES, AND LUNGS. AVOID FORMATION AND BREATHING OF AEROSOLS OR VAPORS.

LABORATORY OPERATIONS SHOULD BE CONDUCTED IN A FUME HOOD, GLOVE BOX, OR VENTILATED CABINET.

AVOID SKIN CONTACT: IF EXPOSED, WASH WITH SOAP AND WATER. AVOID WASHING WITH SOLVENTS. AVOID RUBBING OF SKIN OR INCREASING ITS TEMPERATURE.

BENZENE IS FLAMMABLE AND EXPLOSIVE. KEEP AWAY FROM SPARKS AND OPEN FLAMES. IN CASE OF FIRE, USE CARBON DIOXIDE OR DRY CHEMICAL EXTINGUISHER.

FOR EYE EXPOSURE, IRRIGATE IMMEDIATELY WITH LARGE AMOUNTS OF WATER. FOR INGESTION, DO NOT INDUCE VOMITING. DRINK MILK OR WATER. REFER FOR GASTRIC LAVAGE. FOR INHALATION, REMOVE VICTIM PROMPTLY TO CLEAN AIR. ADMINISTER RESCUE BREATHING IF NECESSARY. REFER TO PHYSICIAN.

IN CASE OF LABORATORY SPILL, WEAR PROTECTIVE CLOTHING DURING CLEANUP. AVOID SKIN CONTACT OR BREATHING OF AEROSOLS OR VAPORS. USE ABSORBENT PAPER TO MOP UP SPILL. AFTER THE RESIDUE HAS EVAPORATED, WASH DOWN AREA WITH WATER. DISPOSE OF WASTE SOLUTIONS AND MATERIALS APPROPRIATELY. MONITOR LABORATORY AIR AND CHECK FOR BENZENE RESIDUES AFTER CLEANUP.

Issued: 10/88

Prepared by the Environmental
Control and Research Program

A. Background

Benzene is a clear, colorless, volatile, highly flammable liquid. It is moderately toxic in animals and man and is carcinogenic, embryotoxic, and mutagenic. Benzene is used in the manufacture of many organic compounds, artificial leather, varnishes, lacquer, and other products. However, because of stringent governmental restrictions on occupational exposure to benzene, its industrial use is being curtailed; where possible, substitute solvents are sought.

The permissible exposure limit is 1 ppm as an 8-hour time-weighted average; the short-term exposure limit is 5 ppm averaged over any 15 minute period (Federal Register, 1987). The current threshold limit value for benzene is 10 ppm (ACGIH, 1988). IARC (1982), p. 97 lists official exposure limits imposed by various countries.

Recent reviews of the chemical and biological properties of benzene include IARC (1982), Mehlman (1983), and Aksoy (1988a).

B. Chemical and Physical Data

1. Chemical Abstract No.: 71-43-2

2. Synonyms:^A

Phene

Pyrobenzole

Benzol

Phenylhydride

Benzole

Cyclohexatriene

Pyrobenzol

Bicarburet of hydrogen

3. Molecular


formula:

C_6H_6

structure:

weight:

78.11

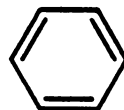

4. Density: Liquid, 0.8787 g/cm^3 at 15°C relative to water at 4°C ; vapor, 2.7 (air = 1.0).

5. Absorption spectroscopy: UV (C_2H_5OH): λ ($\log \epsilon$) = 243 (2.2), 249 (2.3), 256 (2.4), and 261 (2.2). Infrared (Robinson, 1974) and mass (Monahan and Stanton, 1962) spectra have been reported.

"Benzin," "benzine," and "petroleum benzin" are terms used for a mixture of low-boiling aliphatic hydrocarbons and should not be confused with benzene.

6. Volatility: Vapor pressure, 74.6 mm Hg at 20°C. (For vapor pressures below and above 1 atmosphere, see pages D-201 and D-212, respectively, in Weast, 1988.)
7. Solubility: Water, 0.8 parts per 1,000 by weight at 20°C; miscible with ethanol, ether, acetone, chloroform, carbon disulfide, glacial acetic acid, and oils.
8. Description: Clear, colorless, volatile liquid.
9. Boiling point: 80.1°C.
Melting point: 5.5°C.
10. Stability: Stable at normal temperatures; highly flammable.
11. Chemical reactivity: Benzene has relatively low reactivity; it is reduced (to cyclohexane) and oxidized (ring cleavage) only under very stringent conditions. It is, however, subject to substitution reactions such as nitration, sulfonation, halogenation, alkylation, and acylation in the presence of suitable catalysts. For a general review, see standard textbooks (e.g., p. 337 ff, Morrison and Boyd, 1973).
12. Flash point: -11.1°C (closed cup).
13. Autoignition temperature: 538°C.
14. Explosive limits in air: 1.4 - 7.1%.

Fire, Explosion, and Reactivity Hazard Data

1. Use dry chemical fire extinguisher. Fire-fighting personnel should wear air-supplied respirators will full-face masks.
2. Benzene is flammable and its vapors in air can form explosive mixtures (see above).
3. No other conditions contributing to instability are known.
4. No hazardous decomposition products have been identified.
5. Do not expose to sparks or open flames. Use nonspark tools and equipment.

Operational Procedures

The NIH Guidelines for the Laboratory Use of Chemical Carcinogens describe operational practices to be followed when potentially carcinogenic chemicals are used in NIH laboratories. The NIH Guidelines should be consulted to identify the proper use conditions required and specific controls to be implemented during normal and complex operations or manipulations involving benzene.

It should be emphasized that this data sheet and the NIH Guidelines are intended as starting points for the implementation of good laboratory practices when using this compound. The practices and procedures described in the following sections pertain to the National Institutes of Health and may not be universally applicable to other institutions. Administrators and/or researchers at other institutions should modify the following items as needed to reflect their individual management system and current occupational and environmental regulations.

Benzene penetrates various glove materials (Sansone and Tewari, 1978). This factor should be taken into account when handling materials or equipment containing benzene.

1. Chemical inactivation: No validated method reported.
2. Decontamination: Turn off equipment that could be affected by benzene or the materials used for cleanup. If there is any uncertainty regarding the procedures to be followed for decontamination, call the NIH Fire Department (dial 116) for assistance. Use absorbent paper to mop up spill. After the residue has evaporated, wash with copious quantities of water. Glassware should be rinsed (in a hood) with ethanol, followed by soap and water. Animal cages should be washed with water.
3. Disposal: No waste streams containing benzene shall be disposed of in sinks or general refuse. Surplus benzene or chemical waste streams contaminated with benzene shall be handled as hazardous chemical waste and disposed of in accordance with the NIH chemical waste disposal system. Nonchemical waste (e.g., animal carcasses and bedding) containing benzene shall be handled and packaged for incineration in accordance with the NIH medical-pathological waste disposal system. Potentially infectious waste (e.g., tissue cultures) containing benzene shall be packaged for incineration, as above. Burnable waste (e.g., absorbent bench top liners) minimally contaminated with benzene shall be handled as potentially infectious waste and packaged for incineration, as above. Absorbent materials (e.g., associated with spill cleanup) grossly contaminated shall be handled in accordance with the chemical waste disposal system. Radioactive waste containing benzene shall be handled in accordance with the NIH radioactive waste disposal system.
4. Storage: Store in sealed ampoules or in bottles with caps with polyethylene cone liners inside a sealed secondary container. This should be kept in a solvent storage cabinet or explosion-safe refrigerator. Store working quantities of benzene and its solutions in an explosion-safe refrigerator in the work area.

E. Monitoring and Measurement Procedures Including Direct Field Measurements and Sampling for Subsequent Laboratory Analysis

1. Sampling: The procedure for air sampling recommended by the Intersociety Committee on Methods for Air Sampling and Analysis (Sawicki et al., 1975) and NIOSH (1977) is adsorption on

charcoal and desorption with carbon disulfide. Because of the toxicity and flammability of CS₂, as well as its unsuitability when electron capture or photoionization detectors are used in conjunction with gas chromatography, silica gel or commercial polymers with elution by methanol or ethanol have been proposed as substitutes (van Tassel et al., 1981). Methods for the extraction of blood (Snyder et al., 1975) and tissue (Snyder et al., 1977) for subsequent gas chromatographic analysis have been published.

2. **Analysis:** Methods of analysis have been reviewed (Snyder, 1977; Ikeda and Okuno, 1988); IARC (1982, p. 107) provides an overview of methods through 1980. The method of choice is gas chromatography with a variety of detectors. Special high-pressure liquid chromatographic methods have been developed for the separation and quantitation of benzene metabolites for pharmacokinetic studies (Greenlee et al., 1981; Sabourin et al., 1988a).

Biological Effects (Animal and Human)

1. **Absorption:** Benzene is absorbed readily by inhalation, ingestion, and parenteral injection. Percutaneous absorption used to be considered as slight but has recently been shown (in hairless mice) to contribute a significant portion of the total absorbed dose (Susten et al., 1985).
2. **Distribution and pharmacokinetics:** On inhalation of benzene by animals, the 40-50% amount retained (the rest being exhaled unchanged) is transported by the blood (preferentially, the red cells) to the liver, bone marrow, and tissues with high lipid content where it is subsequently metabolized (Schrenk et al., 1941). Excretion in the lungs follows a biphasic pattern, suggesting a two-compartment model for distribution, in the rat (Rickert et al., 1979). Similar results were found in humans where, however, a three-compartment model is indicated (Srbova et al., 1950; Berlin et al., 1979).
3. **Metabolism and excretion:** Schemes of the metabolic pathways of benzene have been depicted (Pellack-Walker et al., 1985; Kalf et al., 1987; and Cooper and Snyder, 1988 among others). Essentially they consist of oxidation to phenol (probably via benzene epoxide) and further by two pathways either to hydroquinone and benzoquinone, or to catechol and 1,2,4benzene triol. These are excreted in the urine as glucuronides and sulfates. In addition, benzene is also excreted in the urine unchanged, the relative amounts of benzene and metabolites depending on the amount administered since the capacity of mammalian liver to oxidize benzene is limited (Sabourin et al., 1987, 1988b). A further excretion product, transmuconic acid, is the result of oxidative opening of the catechol ring. This oxidation has been demonstrated in liver and bone marrow and is mediated by cytochrome P-450 or reactive oxygen species.

Further details are discussed in IARC (1982) and the above-named references.

4. Toxic effects: The acute oral LD50 of benzene in the mouse and rat has been reported as 3-5 g/kg. The inhalation LC50 in these species is 10,000 ppm for a 7-hour exposure. The target organ for acute poisoning is mainly the central and peripheral nervous system, and this poisoning results in muscle tremors, convulsions, salivation, nystagmus, and asphyxiation due to paralysis of the respiratory center. Chronic poisoning effects in man have been summarized as follows: "The net effect in the human is a disease characterized by gradual decreases in white cell levels, among which the granulocytes are most severely affected, but decreases in lymphocytes, the form more sensitive in rodents, are also affected. Red cell and platelet levels may also fall ... Reason for pancytopenia is severe bone marrow damage and the disease is called aplastic anemia" (Snyder, 1984).

These toxic effects have been reviewed in detail by van Raalte and Grasso (1982) and Aksoy (1988b). The biochemical basis for these toxic actions appears to be covalent binding of benzene metabolites to nucleic acids and protein in liver, and inhibition of thymidine uptake in bone marrow DNA (Cooper and Snyder, 1988).

5. Carcinogenic effects: IARC (1982) summarizes its extensive review of the literature as follows: "There is limited evidence that benzene is carcinogenic in experimental animals... There is sufficient evidence that benzene is carcinogenic to man." The evidence has been further reviewed and updated by Aksoy (1988c). Oral administration or inhalation of benzene in BD Wistar rats and Swiss mice results in carcinomas of the Zymbal glands, oral and nasal cavity, skin, forestomach and mammary glands, hepatomas and pulmonary tumors (Maltoni et al., 1985).
6. Mutagenic and teratogenic effects: Benzene is not mutagenic in the Ames test or in the E. coli DNA repair test (DeFlora et al., 1984; Rexroat and Probst, 1985). Evidence to date indicates that benzene is not teratogenic in doses which are fetotoxic and embryo-lethal. In vitro studies on human lymphocytes have demonstrated sister chromatid exchanges on exposure to benzene, but potency in this respect is far lower than that of benzene metabolites (Erexson et al., 1985).

Emergency Treatment

1. Skin and eye exposure: For skin exposure, remove contaminated clothing and wash skin with soap and water. Skin should not be rinsed with organic solvents. Since benzene is readily absorbed through the skin, avoid rubbing of skin or increasing its temperature. For eye exposure, irrigate immediately with copious quantities of running water for at least 15 minutes. Obtain ophthalmological evaluation.

2. Ingestion: Drink plenty of water or milk. Do not induce vomiting. (Vomiting might re-expose the mouth and esophagus.) Refer for gastric lavage.
3. Inhalation: Remove victim promptly to clean air. Administer rescue breathing if necessary.
4. Refer to physician at once. Consider treatment for pulmonary irritation.

References

- ACGIH. 1988. Threshold Limit Values and Biological Exposure Indices for 1988-1989. American Conference of Governmental Industrial Hygienists, Cincinnati, OH.
- Aksoy, M. 1988a. Benzene Carcinogenicity, CRC Press, Inc., Boca Raton, FL.
- Aksoy, M. 1988b. Benzene Hematotoxicity. Ch. 5 in Aksoy, M., 1988a, loc cit.
- Aksoy, M. 1988c. Benzene Carcinogenicity. Ch. 6 in Aksoy, M., 1988a, loc cit.
- Berlin, M., S. Holm, P. Knutsson, and A. Tunek. 1979. Biological threshold limits for benzene based on pharmacokinetics of inhaled benzene in man. Arch Toxicol, Suppl 2:305-310.
- Cooper, K.R. and R. Snyder. 1988. Metabolism. Ch. 4 in Aksoy, M., 1988a, loc cit.
- DeFlora, S., P. Zanicchi, A. Camoirano, C. Bennicelli, and G.S. Badolati. 1984. Genotoxic activity and potency of 135 compounds in the Ames reversion test and in a bacterial DNA-repair test. Mutat Res 133:161-198.
- Erexson, G.L., J.L. Wilmer, and A.D. Kligerman. 1985. Sister chromatid exchange induction in human lymphocytes exposed to benzene and its metabolites in vitro. Cancer Res 45:2471-2477.
- Federal Register. 1987. OSHA Occupational Exposure to Benzene; Final Rule. 52(No. 176, 11 Sept):34562-34578.
- Greenlee, W.F., E.A. Gross, and R.D. Irons. 1981. Relationship between benzene toxicity and the disposition of ¹⁴C-labeled benzene metabolites in the rat. Chem Biol Interact 33:285-299.
- IARC, International Agency for Research on Cancer. 1982. IARC Monographs on the Evaluation of the Carcinogenic Risk of Chemicals to Humans. 29:93-148, Lyon, France.
- Ikeda, M. and T. Okuno. 1988. Analytical techniques. Ch. 2 in Aksoy, M., 1988a, loc cit.
- Kalf, G.F., G.B. Post, and R. Snyder. 1987. Solvent toxicology: Recent advances in the toxicology of benzene, the glycol ethers, and carbon tetrachloride. Ann Rev Pharmacol Toxicol 127:399-427.
- Maltoni, C., B. Conti, G. Cotti, and F. Belpoggi. 1985. Experimental studies on benzene carcinogenicity at the Bologna Institute of Oncology: Current results and ongoing research. Am J Ind Med 7:415-446.
- Mehlman, M.A. (ed). 1983. Carcinogenicity and Toxicity of Benzene. Adv Med Env Toxicol, Vol. 4. Princeton Scientific Publishers, Inc., Princeton, NJ.

- Monahan, J.E. and H.E. Stanton. 1962. Mass spectra resulting from high-energy electron impact on some hydrocarbon molecules. *J Chem Phys* 37:2654-2661.
- Morrison, R.T. and R.N. Boyd. 1973. *Organic Chemistry*, 3rd ed. Allyn and Bacon, Inc., Boston, MA.
- NIOSH. 1977. *NIOSH Manual of Analytical Methods*, 2nd ed., Vol. 3. Method S311-1-S311-8.
- Pellack-Walker, P., J.K. Walker, H.H. Evans, and J.L. Blumer. 1985. Relationship between the oxidation potential of benzene metabolites and their inhibitory effect on DNA synthesis in L5178YS cells. *Med Pharmacol* 28:560-566.
- Rexroat, M.A. and G.S. Probst. 1985. Mutation tests with *Salmonella* using the plate-incorporation assay. *Prog Mutat Res* 5:201-212.
- Rickert, D.E., T.S. Baker, J.S. Bus, C.S. Barrow, and R.D. Irons. 1979. Benzene disposition in the rat after exposure by inhalation. *Toxicol Appl Pharmacol* 49:417-423.
- Robinson, J.W. (ed). 1974. *CRC Handbook of Spectroscopy*. Vol. 2, page 25, CRC Press, Cleveland, OH.
- Sabourin, P.J., B.T. Chen, G. Lucier, L.S. Birnbaum, E. Fisher, and R.F. Henderson. 1987. Effect of dose on the absorption and excretion of [^{14}C]benzene administered orally or by inhalation in rats and mice. *Toxicol Appl Pharmacol* 87:325-336.
- Sabourin, P.J., W.E. Bechtold, and R.F. Henderson. 1988a. A high pressure liquid chromatographic method for the separation and quantitation of water-soluble radiolabeled benzene metabolites. *Anal Biochem* 170:316-327.
- Sabourin, P.J., W.E. Bechtold, L.S. Birnbaum, G. Lucier, and R.F. Henderson. 1988b. Differences in the metabolism and disposition of inhaled [^3H]benzene by F344/N rats and B6C3F mice. *Toxicol Appl Pharmacol* 94:128-140.
- Sansone, E.B. and Y.B. Tewari. 1978. The permeability of laboratory gloves to selected solvents. *Am Ind Hyg Assoc J* 39:169-174.
- Sawicki, E. (Chairman). 1975. Organic solvent vapor in air: Analytical method. *Health Lab Sci* 12:394-402.
- Schrenk, H.H., W.P. Yant, S.J. Pearce, F.A. Patty, and R.R. Sayers. 1941. Absorption, distribution and elimination of benzene by body tissues and fluids of dogs exposed to benzene vapor. *J Ind Hyg Toxicol* 23:20-34.
- Snyder, C.A. 1977. Analytical techniques. In Laskin, S., and B.D. Goldstein. Benzene toxicity, a critical evaluation. *J Toxicol Environ Health* [Suppl] 2:5-22.
- Snyder, C.A., M.N. Erlichman, B.D. Goldstein, and S. Laskin. 1977. An extraction method for the determination of benzene in tissue by gas chromatography. *Am Ind Hyg Assoc J* 38:272-276.
- Snyder, C.A., S. Laskin, and B.D. Goldstein. 1975. An extractive method for the determination of benzene in blood by gas chromatography. *Am Ind Hyg Assoc J* 36:833-836.
- Srbova, J., J. Teisinger, and S. Skramovsky. 1950. Absorption and elimination of inhaled benzene in man. *Arch Ind Hyg Occup Med* 2:1-8.

- Susten, A.S., B.L. Dames, J.R. Burg, and R.W. Niemeier. 1985. Percutaneous penetration of benzene in hairless mice: An estimate of dermal absorption during tire-building operations. Am J Ind Med 7:323-335.
- van Raalte, H.G.S. and P. Grasso. 1982. Hematological, myelotoxic, clastogenic, carcinogenic, and leukemogenic effects of benzene. Regul Toxicol Pharmacol 2:153-176.
- van Tassel, S., N. Amalfitano, and R.S. Narang. 1981. Determination of arenes and volatile haloorganic compounds in air at microgram per cubic meter levels by gas chromatography. Anal Chem 53:2130-2135.
- Weast, R.C. (ed). 1988. Handbook of Chemistry and Physics, 69th ed. CRC Press, Inc., Boca Raton, FL.